

**2008 - Anaesthesia Circuit Paediatric**

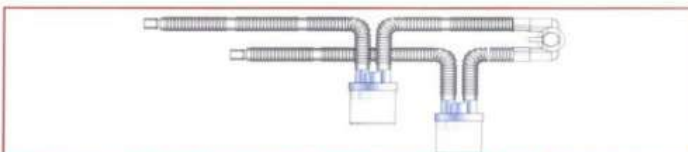
Set comprises : 10mm Corrugated Plastic Hose X 1.6 mtrs (2037), 15mm/22mm Paedi Angle Mount W/Swivel (2062), 15mm Female Hose Mount for 10mm Hose (9065F), 1/2Ltr. Breathing Bag (8201), 15mm/15mm Male/Male Circuit Adaptor (9041), Bag Bleed Valve (8221), 15mm Corrugated Plastic Hose X 0.4 mtrs (2038)

**2019 - Anaesthesia Circuit Adult With 2 Port Wye**

Set comprises : 22mm Corrugated Plastic Hose X 4 mtrs (5037), 2 Port Wye - Adult (2060), 22mm/22mm Male/Female Circuit Adaptors X 2 (9034), 15mm/22mm Angle Mount (9004), 22mm/22mm Male/Male Circuit Adaptor (9033), 2Ltr. Breathing Bag (8204), 22mm Hose Mount - Male (9020M)

**2009 - Ventilator Circuit Plain Adult**

Set comprises : 22mm Corrugated Plastic Hose X 3.2 mtrs (5037), Fixed Wye Adult (9032), 22mm/22mm Male/Female Circuit Adaptors X 2 (9034)

**2012DW - Ventilator Circuit Neonatal With 2 Water Trap**

Set comprises : 10mm Corrugated Plastic Hose X 3.2 mtrs (2037), 1 Port Neonatal Wye (2064), 15mm Female Hose Mount to fit 10mm Hose (9065F) x 2, Paediatric Water Trap (2706) x 2.

Magill Semi - Open Circuits**2013 - (APD - 22) - 22mm All Plastic Disposable**

Set comprises : 22mm Corrugated Plastic Hose X 1.6 mtrs (5037), 22mm/22mm Male/Male Circuit Adaptor (9033), 22mm/22mm Male/Female Circuit Adaptor (9034), 22mm Expiratory Valve (9200), 22mm/22mm Male/Female Angle Mount (9003), 22mm Bag Mount 'T' Piece (9023), 22mm Male Hose Mount (9020M), 2 Ltr. Breathing Bag (8204)

2003 - (PRR - 22) - 22mm Plastic & Rubber Reusable

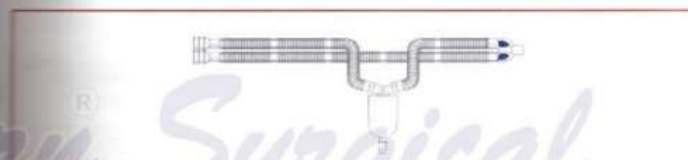
Set comprises : 22mm Corrugated Rubber Hose (8107), 22mm Male/Female Hose Mount Pair (9020), 22mm Expiratory Valve (9200), 22mm/22mm Male/Female Angle Mount (9003), 22mm Bag Mount 'T' Piece (9023), 22mm Male Hose Mount (9020M), 2 Ltr. Breathing Bag (8204)

2006 - (PRR - 23) - 23mm Plastic & Rubber Reusable

Set comprises : 22mm X 105cm Corrugated Rubber Hose (8107), 23mm Male/Female Hose Mount Pair (9019), 23mm Expiratory Valve (9201), 22mm/23mm Male/Female Angle Mount (9001), 23mm Bag Mount 'T' Piece (9025), 22mm Male Hose Mount (9020M), 2 Ltr. Breathing Bag (8204)

2026 - (MRR - 23) - 23mm Metal & Rubber Reusable

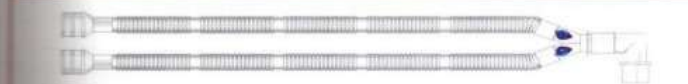
Set comprises : 22mm X 105cm Corrugated Rubber Hose (8107), 23mm Expiratory Valve BC (9202), 22mm/23mm Male/Male Angle Mount BC (9002), 23mm Bag Mount 'T' Piece W/Lever BC (9027), 2 Ltr. Breathing Bag (8204)

**2009W - Ventilator Circuit Adult With 1 Water Trap**

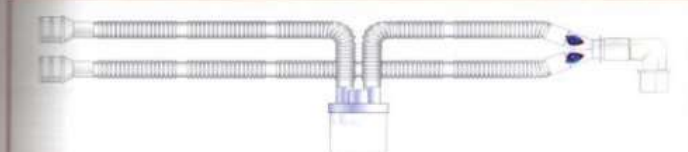
Set comprises : 22mm Corrugated Plastic Hose X 3.2 mtrs (5037), 2 Port Wye - Adult (2060), 22mm/22mm Male/Female Circuit Adaptors X 2 (9034), Water Trap (2701)

2009DW - Ventilator Circuit Adult With 2 Water Trap & Limb

Set comprises : 22mm Corrugated Plastic Hose X 3.2 mtrs (5037), 2 Port Wye - Adult (2060), 22mm/22mm Male/Female Circuit Adaptors X 2 (9034), Water Trap (2701) x 2

**2011 - Ventilator Circuit - Paediatric**

Set comprises : 15mm Corrugated Plastic Hose X 3.2 mtrs (2038), 2 Port Paedi Wye (2063), 15mm/22mm Angle Mount (9004), 22mm Female Catheter Mount X 2 (9008)

**2011W - Ventilator Circuit - Paediatric with 1 water trap**

Set comprises : 15mm Corrugated Plastic Hose X 3.2 mtrs (2038), 2 Port Paedi Wye (2063), 15mm/22mm Angle Mount (9004), 22mm Female Catheter Mount X 2 (9008) & Paediatric Water Trap (2706)

2011DW - Ventilator Circuit - Paediatric with 2 water trap

Set comprises : 15mm Corrugated Plastic Hose X 3.2 mtrs (2038), 2 Port Paedi Wye (2063), 15mm/22mm Angle Mount (9004), 22mm Female Catheter Mount X 2 (9008) & Paediatric Water Trap (2706) x 2.

Magill Closed Circuits**2014 - (APD - 22) - 22mm All Plastic Disposable**

Set comprises : 22mm Corrugated Plastic Hose X 3.2 mtrs (5037), 22mm/22mm Male/Female Circuit Adaptor X 2 (9034), Dr. Nawathe's Swivel with expiratory valve (9050)

2004 - (PRR - 22) - 22mm Plastic & Rubber Reusable

Set comprises : 22mm X 105cm Corrugated Rubber Hose X 2 (8107), 22mm Male Hose Mount X 2 (9020M), Dr. Nawathe's Swivel with expiratory valve (9050)

2007 - (PRR - 23) - 23mm Plastic & Rubber Reusable

Set comprises : 22mm X 105cm Corrugated Rubber Hose X 2 (8107), 23mm Male Hose Mount X 2 (9019M), Dr. Nawathe's Swivel with expiratory valve (9050)

2027 - (MRR - 23) - 23mm Metal & Rubber Reusable

Set comprises : 22mm X 105cm Corrugated Rubber Hose X 2 (8107), 23mm Male Angled Hose Mount BC X 2 (9047M), 23mm Metal Swivel with expiratory valve (9033)

**3030 - Dr. P. K. Sen's Weaner (Recovery T) with tubing**

Used for rapid oxygenation of patients' blood after surgical procedures. Will shorten the time spent per patient in the recovery room. Fits all standard endotracheal tubes and tracheostomy tubes (15mm). Allows supplemental oxygen therapy in the recovery room. Oxygen is delivered through a small bore tubing connected to the oxygen barb on the weaner. A 6" piece of 22mm tubing acts as an Oxygen Reservoir. A Four-way Tee with Cap and Loop allows suction procedures to be accomplished without removing the weaner.

* Carrying Warranty at Rajkot Office only on site warranty not provided free of cost.

* Due to constant upgradation design, price and features are subject to change any time without prior notice.



Western Surgical

WESTERN SURGICAL
EQUIPMENTS PRIVATE LIMITED

25

* www.westernsurgical.in *

Email : westernsurgical@gmail.com

Ph. +91-281-2232871 / 2243159

Mo. +91 94280 03750/52/54

© Copyright Reserved by : Western Surgical - Rajkot.